

# COMPANION 5™

5 LPM At-Home Oxygen Concentrator  
CE Certified for International Distribution



## Quality. Performance. Reliability.

The Companion 5's small profile fits easily into any home and stays out of the way of your everyday activities. Built with robust internal components, the Companion 5 also has a unique air intake system for maximum durability while producing up to 5 LPM. The proven quality performance maximizes patient confidence and minimizes maintenance calls.

And now, thanks to the optional telehealth solution—myCAIRE™—your provider can periodically monitor this device's location and performance, as well as troubleshoot alerts—all from the convenience of a desktop computer or smart device.



### QUALITY TECHNOLOGY

autoFLOW™ Technology keeps your power consumption low, saves you money, and pays for itself in energy savings.



### PEACE OF MIND PERFORMANCE

Designed with an alert system and easy-to-understand controls—the Companion 5 comes equipped with wireless technology which can connect you to myCAIRE—a telehealth solution that your provider can use to help ensure you get the support you need.



### COMFORT AND EASE OF USE

Stream your favorite TV show, cheer on your sports team, or spend time with loved ones, and never worry about your oxygen delivery. The Companion 5 is perfect to help you go about your daily routine when you are at home.



#### CAIRE Contact Information:

Providers Toll Free

1.800.482.2473

Patients Toll Free

1.877.704.0878

[www.caireinc.com](http://www.caireinc.com)

# COMPANION 5™

5 LPM At-Home Oxygen Concentrator  
CE Certified for International Distribution

## Includes

- Oxygen Concentrator
- Humidifier Adaptor Tube

## Product Features

- Small profile fits easily into home
- Low power draw
- Easy-access humidifier stand
- Simple-to-read controls and alarms
- CE certified
- Optional telehealth solution — myCAIRE



## Specifications

**FLOW RATE**  
0.5 - 5.0 LPM

**SOUND LEVEL**  
43.74 dB(A) ±0.45 dB(A)  
at 3 LPM\*

**MAXIMUM TUBING**  
50 ft (15.2 m)

**DIMENSIONS**  
21.5 x 12.5 x 13.5 inches  
(54.6 x 31.8 x 34.2 cm)

**OXYGEN CONCENTRATION**  
90% +5.5%/-3%

**WARRANTY**  
3 Years

**WEIGHT**  
36 lbs (16.3 kg)

**POWER CONSUMPTION**  
285 W @ 2 LPM,  
350 W Maximum



## Available Accessories



### HUMIDIFIER ADAPTOR & BOTTLE

Humidifier Adaptor (9 in/23 cm)


Ask your doctor or therapist if humidification is right for you.



### NASAL CANNULAS

Available in 7 ft (2.1 m) and 25 ft (7.6m)

Talk to our Customer Service team to learn more.

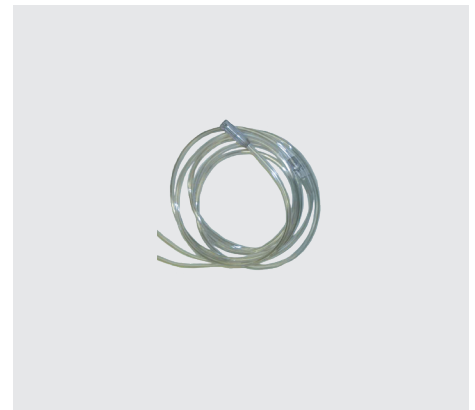

### OXYGEN TUBING

Comes in: 7 ft (2.1 m), 25 ft (7.6 m), and 50 ft (15.2 m).



**CAIRE Contact Information:**  
Providers Toll Free  
1.800.482.2473  
Patients Toll Free  
1.877.704.0878  
[www.caireinc.com](http://www.caireinc.com)

\*per test standard Nr. 14-110/2018 MDS-Hi  
SEE PRODUCT WARRANTY STATEMENT FOR COMPLETE INFORMATION. Please consult the applicable product instructions for use for product indications, contraindications, warnings, precautions, and detailed safety information. myCAIRE is a telehealth solution intended to share specified device data from the enabled CAIRE oxygen concentrator. myCAIRE does not affect oxygen concentrator settings, performance, or oxygen therapy, and is not intended for patient monitoring, diagnosis or treatment. Please seek a doctor's advice before making any medical decisions.  
© Copyright 2021 CAIRE Inc. All Rights Reserved.

ML-CONC0173 D